



Accurate SHBG Results with Test Cup Convenience

*Quality Results
for Quality Life*




Using Total Testosterone and SHBG, results can be provided as Free Androgen Index (FAI), Free testosterone*, Weakly bound testosterone* and Bioavailable testosterone*.

**Require Albumin measurement.*

TOSOH BIOSCIENCE

Available for Labs and Workloads of All Sizes

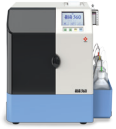

AIA-360
36 Tests Per Hour



AIA-900 Benchtop
90 Tests Per Hour



AIA-900 Loader Model
90 Tests Per Hour




AIA-900 with 9 tray sorter
90 Tests Per Hour
Also available with 19 tray sorter



AIA-2000
200 Tests Per Hour

Unit Dose Test Cup Reagents

- Single, unitized cups use a dry reagent format that requires no pre-mixing, no pre-measuring, no on-board refrigeration and no waste
- Dry reagent format ensures calibration stability of up to 90 days for most assays
- Reagents have a one year shelf life from date of manufacture
- The AIA test cups and trays are labeled with code and lot number for automated scheduling and inventory management on the AIA-2000 and AIA-900 sorter models
- ST AIA-PACK reagents provide a faster turnaround time due to a 10 minute incubation protocol
- Extensive test menu provides a wide range of immunoassay tests including thyroid, cardiac markers, tumor markers, reproductive hormones, anemia markers, metabolic markers, diabetes, kidney markers and more

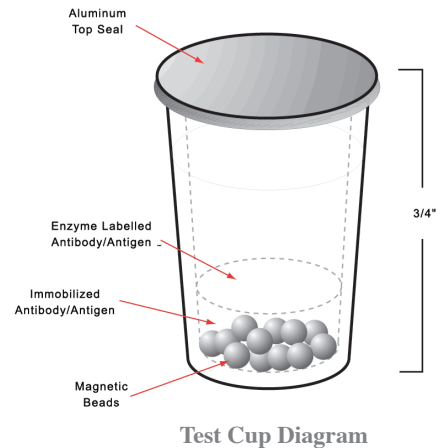

Test Cup Diagram



TOSOH BIOSCIENCE

TOSOH

Tosoh Bioscience, Inc.
6000 Shoreline Court, Suite 101
South San Francisco, CA 94080
Tel: (800) 248-6764
Fax: (650) 615-0415
www.tosohbioscience.us

"AIA-PACK" and "AIA" are the registered trademarks of Tosoh Corporation in United States, Canada, etc.