

Minimize the Risk of False Lab Test Results



Tosoh is Your Simple Solution to Biotin Interference

The increase in biotin supplementation use combined with the limitations in many immunoassays has led to false lab results, misdiagnosis and mismanagement of patients, triggering the FDA to issue a warning against biotin interfering with some lab tests¹⁻³.

Tosoh's ST AIA-PACK[®] Test Menu of immunoassays use a proprietary method that is **free from biotin interference**.


Not all immunoassays are the same. Using the right immunoassay gives you confidence in the result you deliver. **Choose wisely. Choose Tosoh.**

TOSOH BIOSCIENCE

www.tosohbioscience.us



Tosoh's Immunoassays are Free from Biotin Interference

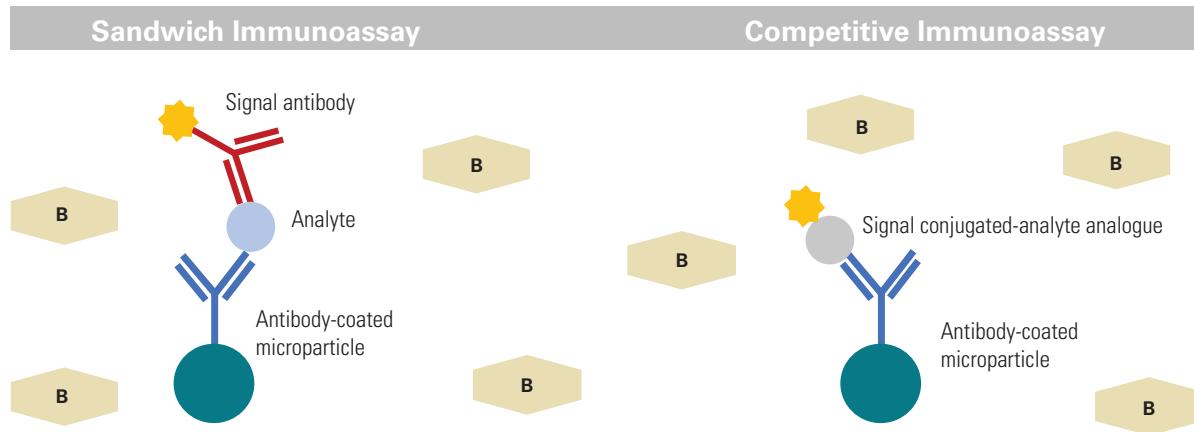


Biotin does not affect Tosoh assays because our assays do not employ the biotin-streptavidin pair to isolate our signal. Our antibodies are coated on microparticles and form the basis of our proprietary dry-reagent **Unit Dose Test Cup Technology**.



Tosoh's Unit Dose Test Cup Technology

Tosoh Bioscience's proprietary dry-reagent, unit dose test cup technology is **free from biotin interference**, offers calibration stability of up to 90 days, and achieves consistency and accuracy with reduced wasted. The test cup is interchangeable across Tosoh's portfolio of automated immunoassay analyzers. Tosoh analyzers are suited for a variety of different throughput requirements, from small physician office labs to large reference labs, and are easy to use, robust, and precise, and offer a simplified workflow with minimal downtime.

Available for Labs and Workloads of All Sizes

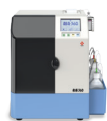

AIA-360



AIA-900 Benchtop



AIA-900 with 9 Tray Sorter
Also available with 19 tray sorter option and as Loader model



AIA-2000

Choose Wisely. Choose Tosoh.

Universal Reagent

Tosoh's AIA-PACK test cup format works with every Tosoh automated immunoassay system allowing for a seamless transition from one system to the other, ensuring consistent results in an efficient and economical process.

Save Time and Money

Tosoh's AIA-PACK test cups are single, unitized cups that use a dry reagent format that ensures calibration stability of up to 90 days.

No Contamination, More Traceability

Because there is no transfer of reagents the risk of contamination is eliminated. The unique bichromatic fluorescence kinetic measurement ensures a high analytical and functional sensitivity for all assays. AIA-PACK test cups and trays are labelled with the assay code and lot number for automated scheduling and inventory.

For more information, call 1.800.248.6764

Tosoh products are for Prescription use only as In-Vitro Diagnostics

TOSOH BIOSCIENCE

www.tosohbioscience.us

"Tosoh Bioscience", "AIA-PACK" and "AIA" are the registered trademarks of Tosoh Corporation
Part # AIA-BIO-SS01 Rev.1/072019

