



# IVF SOLUTIONS



**E2 Results in 18  
to 20 Minutes**

## Why Wait?

*Tosoh Immunoassay Analyzers  
Deliver Accurate Results Quickly*

**Fast** - First Result in 18 to 20 Minutes

**Efficient** - Simple Automated Testing


**Accurate** - Correct Results the First Time

**Consistent** - Common Unit Dose Test Cup Reagent

**Cost Saving** - 90 Day Calibration Stability for all Fertility Assays

**Comprehensive** - Full Immunoassay Menu Available, including a Full Menu of Reproductive Hormone Assays

**Biotin Free** - All Tosoh's AIA-PACK<sup>®</sup> Immunoassays are Free from Biotin Interference

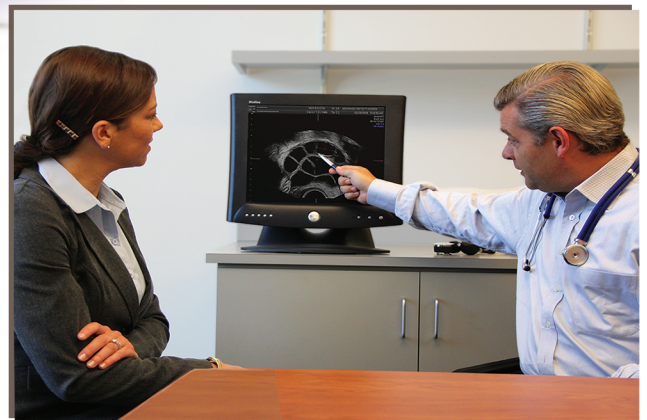

*Tosoh's fast results allow review of ultra-sound  
and E2 on the same visit.*

**TOSOH BIOSCIENCE**

# Unit Dose Test Cup Technology

- Single, unitized cups use a dry reagent format that requires no pre-mixing, no pre-measuring, no on-board refrigeration and no waste
- Dry reagent format ensures calibration stability of up to 90 days for most assays
- Reagents have a one year shelf life from date of manufacture
- The AIA test cups and trays are labeled with code and lot number for automated scheduling and inventory management on AIA-2000 and AIA-900 sorter models
- ST AIA-PACK reagents provide a faster turnaround time due to a 10 minute incubation protocol
- Extensive test menu provides a wide range of immunoassay tests including thyroid, cardiac markers, tumor markers, reproductive hormones, anemia markers, metabolic markers, diabetes, kidney markers and more
- Biotin does not affect Tosoh assays because the methodology does not utilize the biotin-streptavidin pair to isolate the signal



## Available for Labs and Workloads of All Sizes



**AIA-360**  
36 Tests Per Hour




**AIA-900 Benchtop**  
90 Tests Per Hour



**AIA-900 Loader Model**  
90 Tests Per Hour



**AIA-900 with 9 tray sorter**  
90 Tests Per Hour  
*Also available with 19 Tray Sorter*

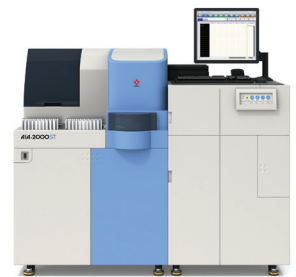

**AIA-2000**  
200 Tests Per Hour

**Choose Wisely. Choose Tosoh.**



**TOSOH BIOSCIENCE**

*"AIA-PACK" and "AIA" are registered trademarks of Tosoh Corporation*

Tosoh Bioscience, Inc.  
6000 Shoreline Court, Suite 101  
South San Francisco, CA 94080  
Tel: (800) 248-6764  
Fax: (650) 615-0415  
[www.tosohbioscience.us](http://www.tosohbioscience.us)